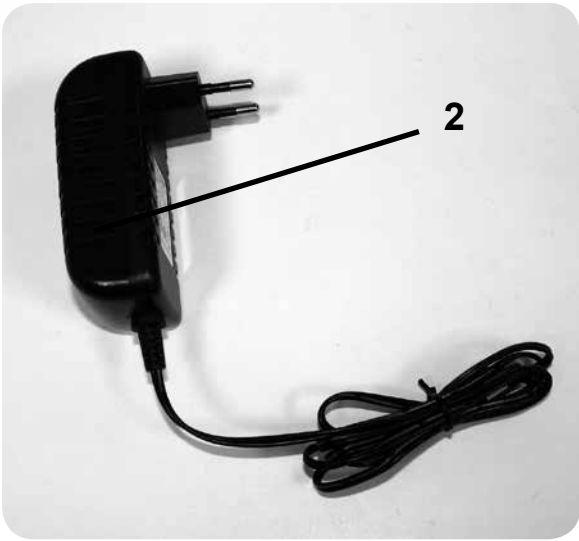
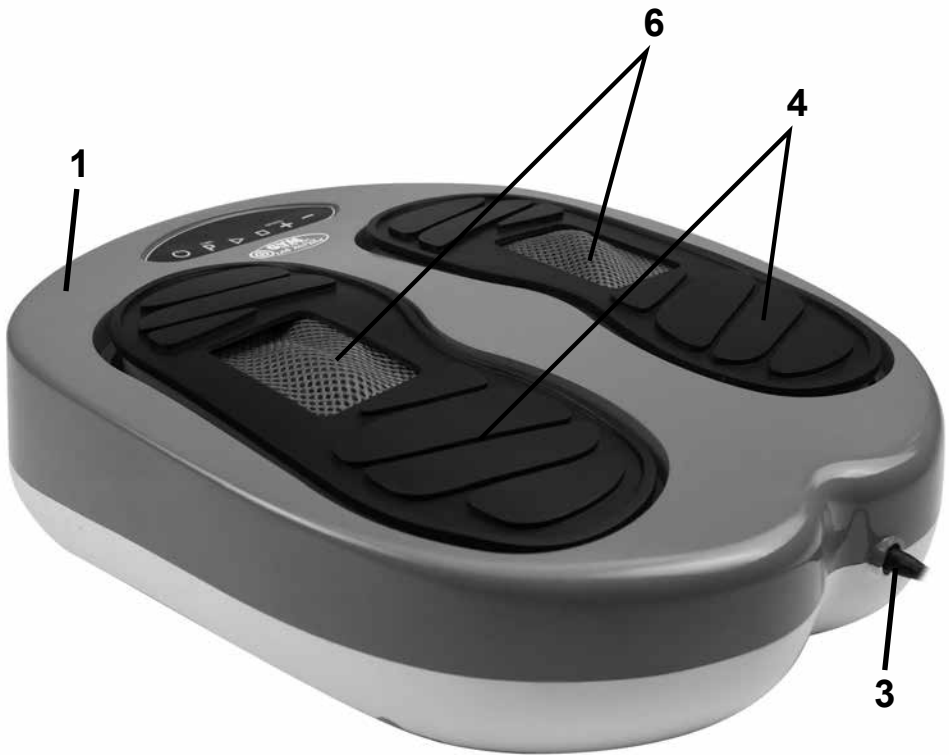
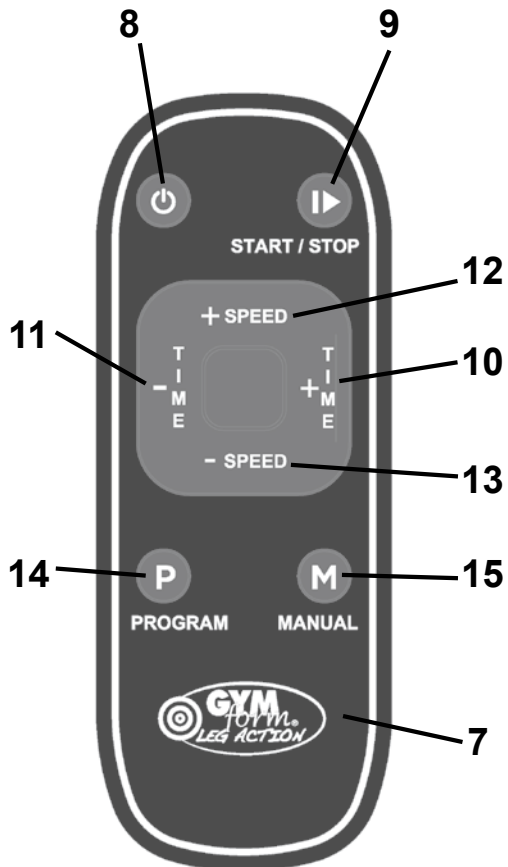
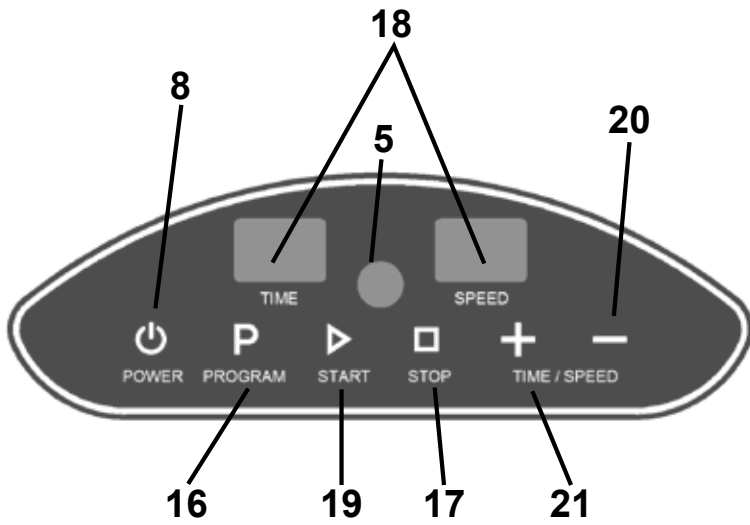
-  *Instruction Manual*
-  *Manual de instrucciones*
-  *Mode d'emploi*
-  *Bedienungsanleitung*
-  *Istruzioni per l'uso*
-  *Manual de instruções*
-  *Gebruiksaanwijzing*

International patents pending. All imitations will be prosecuted.
Gymform® Leg Action is a registered EU/CTM trade mark.

Safety Warnings

Read the instruction manual before using Gymform® Leg Action for the first time.

Always place Gymform® Leg Action on a flat surface.

Ensure that children do not use Gymform® Leg Action without adult supervision.

Do not allow Gymform® Leg Action to make contact with water.

Do not use Gymform® Leg Action if your hands are wet.

Unplug the device after use.

As with any other physical exercise device, consult your doctor before you start using Gymform® Leg Action.

This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.

This product is not intended for medical use. It is intended only for massages.

Consult your physician before using this product if:

You are pregnant

You have a pacemaker, an implant or other aids

You suffer from diabetes

You have a circulatory disorder

You have concerns about your health

The elderly and infirm should consult their doctor before using the product, even if there are no specific pre-existing conditions.

Because skin friction can be painful, it is recommended that the device not be directly applied to the skin. Stop using immediately if you feel discomfort or pain.

Do not use this product at bedtime. The

massager has a stimulating effect and can delay sleep.

To avoid a short circuit, connect the device correctly and keep it away from metal (needles), waste and water.

Use the device on a flat surface. If the device is not stable, it could tip over and cause injury.

Remove the adapter from the power outlet when not in use.

If you suspect a problem, turn off the device as soon as possible.

When not in use for extended periods, store the device in its packaging in a safe, dry, cool place.

Do not place any other objects on the massager.

If the power cord is damaged, it must be replaced, but only by an authorized and qualified person.

To reduce the risk of electric shock, do not use the device on wet skin.

Keep the product away from moisture. Do not use the product if it or its parts are wet or damp.

Keep the device and its cable out of the reach of children under the age of eight.

To avoid injury and irreparable damage, always sit while using the massager. To avoid dangerous injury, never stand on the unit or use it forcefully.

Do not attempt to disassemble or repair the product. Repairs may only be carried out by suitably qualified technicians.

Product Description:

1. Gymform® Leg Action main unit
2. Power adapter
3. Power jack
4. Foot pads
5. Control panel
6. Shiatsu massage rollers
7. Remote control
8. ON/OFF button 
9. Start / Stop button 
10. More Time button (Time+)
11. Less Time button (Time-)
12. More speed button (Speed+)
13. Less speed button (Speed-)
14. Automatic programme button (P) on remote control
15. Manual programme button (M) on remote control
16. Programme button (P) on device
17. Button to stop the device
18. Display panels
19. Button to start the device 

- 20. - Button on device
- 21. + Button on device

The Set Contains:

- 1 Gymform® Leg Action main unit
- 1 Remote control
- 1 Instruction manual
- 1 Exercise guide

Instructions for Use

- Operating features
- Speeds from 1 to 15.
- Six specifically designed programmes.

Note:
The device is designed to be used on the sole of the foot, the legs and back. Do not use it on the waist, neck or other parts of the body not described in this manual. Do not use the device for purposes other than those specified in this manual and the exercise guide.

1. Place the unit (1) on a flat, level surface.
2. Plug one end of the adapter (2) into the mains and the other into the power jack (3).
3. The functions of the Gymform® Leg Action can be selected on the control panel (5) on the device or on the remote control (7), as described below. The remote control runs on two AAA batteries (not included). Insert them into the battery compartment by following the polarity marks.
4. To turn the device on and off, press the ON/OFF button (8) .
5. When the device is on, press the “P” button (14) on the remote control to choose one of the Gymform® Leg Action’s four automatic programmes (P1, P2, P3 and P4). Press the “M” button (15) on the remote control to choose one of the two manual programmes (M0 and M1). To choose the programme on the device, just press the Programme button (P) on the device (16). Press this button to switch from automatic to manual programmes. The chosen setting is shown on the display panels (18).
6. (see Programme Table).

Programme Table:

AUTOMATIC PROGRAMMES

P1 –
Vibration only;
Range of intensity levels: 6-15
Time: 10 minutes

P2 –
Vibration and vibration + massage
1 minute.....vibration
2 minutes.....vibration + massage
3 minutes.....vibration
4 minutes.....vibration + massage
5 minutes.....vibration
6 minutes.....vibration + massage
7 minutes.....vibration
8 minutes.....vibration + massage
9 minutes.....vibration

10 minutes vibration + massage
Range of intensity levels: 9-15
Time: 10 minutes

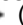
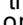
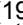
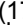
P3 –
Vibration only
Range of intensity levels: 9-15
Time: 10 minutes

P4–
Vibration and vibration + massage
2 minutes.....vibration
2 minutes.....vibration + massage
3 minutes.....vibration
3 minutes.....vibration + massage
Range of intensity levels: 9-15
Time: 10 minutes

MANUAL PROGRAMMES

M0
Manual Programme, vibration only
Range of intensity levels: 1-15
Time: 10 minutes

M1
Manual Programme, vibration + massage
Range of intensity levels: 1-15
Time: 10 minutes

7. To start and stop the chosen programme, press the Start/Stop button  (9).
8. It is not possible to set the time or speed on automatic programmes.
9. If you choose one of the two manual programmes (M0 and M1), you can set a time of 1 to 10 minutes before starting the programme.
 - Press the Less Time button (T-) (11) on the remote control or the - button on the unit (20) to reduce the time.
 - Press the More Time button (T+) (10) on the remote control or the + button on the unit (21) to increase the time.
10. Once you have chosen the time, press the Start/Stop button  (9) on the remote control or the Start button  (19) on the device to start the programme. When the programme begins, you can set the speed (from 1 to 15) at any time.
 - Press the More Speed button (S+) (12) to increase the speed.
 - Press the Less Speed button (S-) (13) to reduce the speed.
11. To turn off the device, first press the Start/Stop button  (9) on the remote control or the stop button on the device  (17). Then press the ON/OFF button on the device to turn it completely off.
12. Disconnect the device from the mains after use.

Positions for using the device:
The device can also be used on the legs and back. Consult the exercise guide on how to use the product.

NOTE: Never apply all your weight on the

device, as this may damage it.

CLEANING

Always unplug the device from the mains before cleaning it.

Wipe the device and foot pads (4) with a damp cloth and mild soap. Do not use abrasive or corrosive cleaners. The foot pads can be removed for easy cleaning (4).

TECHNICAL SPECIFICATIONS

Input: 24,0 V 
Power: 30 W

Adapter Technical Specifications

Input: 100-240 V

Output 24,0 V  1,5 A

Technical information on the FY0422401500E adapter

The product acquired "GYMFORM LEG ACTION, ref. VDPGYCIND0090", is designed to work only with the adapter described in the table(*).

For your own safety and the safety of others, we do not recommend using this adapter with any other device.

(*)

Manufacturer's name or trademark, commercial registration number and address /

PRINCIPAL
PRODUCT: INDUSTEX.
Av. Països Catalans
34, 8th floor. 08950,
Esplugas Llobregat,
Spain.

Adaptador: Shenzhen
Fuyuanjian Power
Co.,Ltd 8FA A9 Building
Tianrui Industrial Park
Fuyong Town baoan
District Shenzhen
Guangdong, China

Model identifier / FY0422401500E

Input voltage / 100-240 V

AC input frequency / 50-60 Hz

Output voltage / 24,0 V

Output current / 1,5 A

Output power / 36,0 W

Average active efficiency / 89,47 %

Efficiency at low load
(10%) / 86,21 %

No-load power
consumption / 0,086 W

AUSTRALIA & NEW ZEALAND ONLY

This product has the benefit of certain Consumer Guarantees. These are prescribed by the Australian Consumer Law & the New Zealand Consumer Guarantees Act 1993 both of which provide protection for consumers. There is no express warranty for this product in Australia or New Zealand. The above paragraph refers to other countries.

WARRANTY

This product is covered by a warranty against manufacturing defects subject to the time periods stipulated by the legislation in force in each country.

ATTENTION DISPOSAL OF MATERIALS

A symbol of a crossed-out wheeled bin means you should find out about and follow local regulations about disposing of this kind of products.

Do not dispose of this product as you would other household waste.

Dispose of this device in accordance with the corresponding local regulations.

Electrical and electronic devices contain hazardous substances that can have harmful effects on the environment and/or human health and should be recycled properly.

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ISL GmbH
Best Direct (International) Benelux B.V.
Best Direct International Ltd,
ISL ITALY S.r.l.